



Ready-To-Run


DNA MARKERS

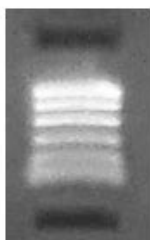


a. 1.2%, 96 wells, 5 min

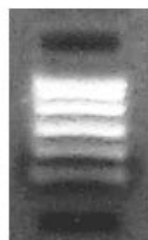

b. 2.2%, 96 wells, 5 min

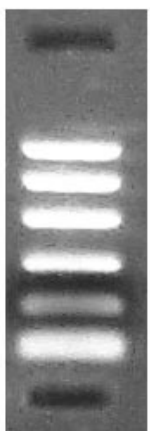

c. 1.2%, 48 wells, 10 min

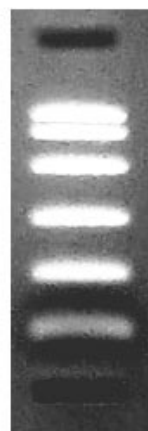

d. 2.2%, 48 wells, 11 min

Fig 1. Ready-To-Run DNA Markers separated on each of the four Ready-To-Run Agarose Gels.

- Prediluted mixture of 6 DNAs is specifically designed for use with Ready-To-Run™ Agarose Gels.
- Supplied at a concentration of 2 µg/ml per DNA band. Marker can be used straight from the vial.
- DNAs of 100, 250, 500, 1 000, 2 000, and 5 000 bp produce easily read patterns on all Ready-To-Run gels.
- Sufficient DNA is supplied for 20 gels at 2 columns of 8 wells or 40 gels at 2 columns of 4 wells.
- Standard buffer contains Bromophenol Blue tracking dye.

Ordering information

product	quantity	code no.
Ready-To-Run DNA Marker	25 µg	80-6485-65
Related products:		
Ready-To-Run Separation Unit	1	80-6460-95
Ready-To-Run Satellite Unit	1	80-6461-14
Ready-To-Run Loading Guide Plates	2	80-6462-09
Ready-To-Run 96-well 1.2% Agarose Gel	10/pk	80-6461-33
Ready-To-Run 48-well 1.2% Agarose Gel	10/pk	80-6461-71
Ready-To-Run 96-well 2.2% Agarose Gel	10/pk	80-6482-80
Ready-To-Run 48-well 2.2% Agarose Gel	10/pk	80-6482-99

Asia Pacific +852 2811 8693, **Australasia** +61 2 9899 0999, **Austria** 01 576 0616 24, **Belgium** 0800 73888, **Canada** 1 800 463 5800, **Central, Eastern & Southeastern Europe** +43 1 982 3826, **Denmark** 4516 2400, **Finland** 358 (0)9 512 3940, **France** 0169 35 67 00, **Germany** 0761 4903 402, **Italy** 02 27322 1, **Japan** 81 3 5331 9336, **Latin America** +55 11 3667 5700, **Middle East and Africa** +30 (1) 96 00 687, **Netherlands** 0165 580 410, **Norway** 23 18 58 00, **Portugal** 21 417 7035, **Russia & other C.I.S. & N.I.S.** +7 (095) 232 0250, **Southeast Asia** +60 3 724 2080, **Spain** 93 594 49 50, **Sweden** 018 612 19 00, **Switzerland** 01 802 81 50, **U.K.** 0800 616 928, **USA** +1 800 526 3593

www.electrophoresis.amershambiosciences.com

Amersham Biosciences UK Limited,
Amersham Place Little Chalfont, Buckinghamshire England HP7 9NA

Amersham Biosciences AB
SE-751 84 Uppsala Sweden

Amersham Biosciences Inc.
800 Centennial Avenue PO Box 1327, Piscataway NJ 08855 USA

Amersham Biosciences Europe GmbH
Munzinger Strasse 9 D-79111 Freiburg Germany

Ready-To-Run is a trademark of Amersham Biosciences Ltd. or its subsidiaries.
Amersham and Amersham Biosciences is a trademark of Amersham plc.

All goods and services are sold subject to the terms and conditions of sale of the company within the Amersham Biosciences group which supplies them. A copy of these terms and conditions is available on request.

© 2001 Amersham Biosciences Inc. All rights reserved.

Printed in the USA

